

Water Distiller






Innova Bio-Meditech is one of the professional solution provider of laboratory and medical devices. Firmly committed to our mission of "sharing innovative bio-meditech solutions with the world", we are dedicated to innovation in the fields of Biology Project, Life Science, Pharmacy Industry and Medical Treatment.

Innova Bio-Meditech possesses a sound distribution and service network with business partners in North and Latin America, Europe, Africa and Asia-Pacific etc. We have built up a well established R&D, manufacture network with 3 centers in Qingdao, Shanghai and Suzhou. Inspired by the needs of our customers, we adopt advanced technologies and transform them into accessible innovation. This means constant effort and research, in order to more fully understand the developments of the market, INNOVA

The passion for science

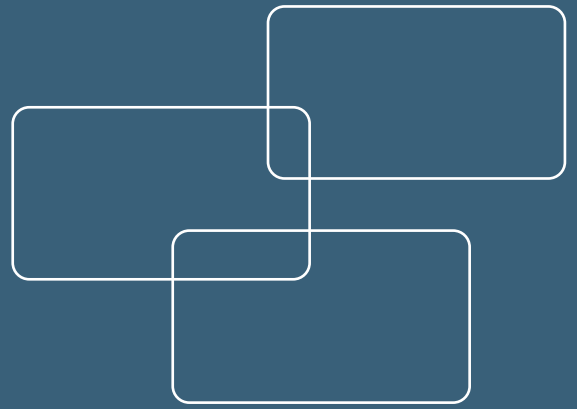

Content

Electric Water Distiller(201)-----	02
Electric Water Distiller(304)-----	03
Automatic cut-off Water distiller(201)-----	04
Double water distiller(304)-----	05
Tower double water distiller(304)-----	07

A scientist wearing a white lab coat, a white hairnet, and blue gloves is looking through a microscope in a laboratory. The scene is lit with a cool blue light, creating a professional and scientific atmosphere. The scientist's hands are positioned to adjust the microscope, and the background shows other laboratory equipment and a blurred figure of another person in a lab coat.

Application

The stainless steel electric automatic distillation Water distiller is especially designed for medicine, plating, lab and sanctuary use. Using tap water, it is safe, stable, capacious, properly structured ,and easy to use.

Electric Water Distiller(201)

Evaporating pot is made of quality stainless steel with all the Argon welding surfaces safe and smooth.



Technical Parameter

Model	IWD-3
Output of water	3L/h
Consumption power	2.5Kw
Mains voltage	220V
Working size	33×24×62(cm)

Electric Water Distiller(304)



IWD-5



IWD-10U



IWD-20U

Features

- All of products are made from stainless steel.
- Condenser is made of seamless stainless-steel pipes and so with very good cooling effects & great capacity of water produce.
- The cooled water returning to the vaporizing drum can save energy and water.
- Drain valve can push the water in container, prevent bilge.
- Immerge electric heating tube made from copper,heat efficiency is better.smooth.

Technical Parameter

Model	IWD-5U	IWD-10U	IWD-20U
Output of water	≥5L/h	≥10L/h	≥20L/h
Consumption power	4.5kW	7.5kW	15kW
Mains voltage	220V/50Hz	380V/50Hz	380V/50Hz
Package size	35×26×78(cm)	42×30×90(cm)	51×36×80(cm)
Gross weight	7kg	10kg	13.5kg
Net weight	5kg	8kg	11.5kg

Automatic cut-off Water distiller(201)

- This product adopts quality stainless steel and is made out of advanced technology and art.
- The cooler is made of seamless stainless steel tube with a high heat-exchange rate and big volume.
- This product adopts imported apparatus, which is highly sensitive and cut off the power when the tank is nearly out of water with a stable performance to protect the user.



Technical Parameter

Model	IAWD-3
Output of water	3L/h
Consumption power	2.5Kw
Mains voltage	220V
Working size	33×24×62(cm)

Automatic cut-off Water distiller(304)



IAWD-5U



IAWD-10U



IAWD-20U

Features

- All of products are made from stainless steel.
- Condenser is made of seamless stainless-steel pipes and so with very good cooling effects & great capacity of water produce.
- The cooled water returning to the vaporizing drum can save energy and water.
- Drain valve can push the water in container, prevent bilge.
- When lack water, will auto stop heating, when the return water reach the normal level, then continue heating.
- Immerge electric heating tube made from copper, heat efficiency is better.

Technical Parameter

Model	IAWD-5U	IAWD-10U	IAWD-20U
Output of water	≥5L/h	≥10L/h	≥20L/h
Consumption power	4.5kW	7.5kW	15kW
Mains voltage	220V/50Hz	380V/50Hz	380V/50Hz
Package size	35×34×78(cm)	42×38×90(cm)	51×36×80(cm)
Gross weight	7kg	10kg	13.5kg
Net weight	5kg	8kg	11.5kg

Double water distiller(304)

Features

- All of products are made from stainless steel.
- With two container, and firm reasonable.
- Condenser is made of seamless stainless-steel pipes and so with very good cooling effects & great capacity of water produce.
- When lack water, will auto stop heating, when the return water reach the normal level, then continue heating.
- Composed of first distilling, second distilling water and heating two control system.
- Immerge electric heating tube made from copper, heat efficiency is better.



Technical Parameter

Model	IDWD-5U	IDWD-10U	IDWD-20U
Output of water	≥5L/h	≥10L/h	≥20L/h
Consumption power	10kW	7.5kW	15kW
Mains voltage	380V/50Hz	380V/50Hz	380V/50Hz
Package size	67×43×86(cm)	42×38×90(cm)	51×36×80(cm)
Gross weight	32kg	38kg	56kg
Net weight	16kg	18kg	34kg

Tower double water distiller(304)

Features

- All of products are made from stainless steel.
- Use high temperature steam to heat, save the energy.
- Use the steam condensation water.
- Seamless stainless-steel pipes. heat efficiency is better.
- Condenser is made of seamless stainless-steel pipes and with good cooling effects & great capacity of water produce.
- Percolation, pushing ammonia, dividing water and steam when distilling, in order to control the quality.



Technical Parameter

Model	ITWD-50U	ITWD-100U	ITWD-200U	ITWD-400U	ITWD-600U
Output of water	≥50L/h	≥100L/h	≥200L/h	≥400L/h	≥600L/h
Pressure of steam supply	0.4Mpa				
Pressure of cooling supply	0.2Mpa				
Steam consumption	≈1:1				
Ratio of water consumption	≤1:14				



Professional Solution Provider of Laboratory/Pharma Equipment

Water Purification System
Freeze Dryer/Lyophilizer
Biosafety Cabinet/Laminar Flow

Biofermentation Solution
Shaking Incubator/Shaker
Cold Storage Solution

Glassware Washer
Autoclave
Climate Chamber



Qingdao Innova Bio-Meditech Co., Ltd.

Add. : No. 176 Jufeng Road, 266199, Qingdao, China

Tel. : +86 532 8789 0634

Email: info@innobiomed.com

web: www.innovabiomed.com